

USER GUIDE

*DermaSeta™*



The Ultimate Body Treatment System





The Ultimate Body Treatment System

NOTES:

---

---

---

---

---

---

---

---

---

---

---

---



TROUBLESHOOTING

SYMPTOMS

PROCEDURES

• BASIC UNIT

Attachments do not rotate:

Make sure the batteries are in the unit and placed in the right way.

Hair removal seems to take longer than at first, or hair is not removed at all:

Feel the surface of the pad; if it feels smooth, then it is time to replace the pad. Follow instructions for the pad replacement.

Attachments come off when the power switch is applied:

Be sure the attachment is installed correctly. When installing, you should hear a soft 'click' that indicates that the pad is secured.

Decreased output is noticed:

Please note that the top performance of **Derma Seta™** can only be ensured with new alkaline batteries.

• DELUXE UNIT

**Derma Seta™** won't turn on:

In order for **Derma Seta™** to turn on, the batteries need to be sufficiently charged. Prior to first use, charge the unit for 4 hours.

Batteries won't charge:

In order to prevent the batteries from complete power drainage, do not allow **Derma Seta™** to remain in an uncharged condition over a large period of time.

THANK YOU

On behalf of the entire team at Thane Direct, we appreciate your recent purchase of **Derma Seta™**.

**Derma Seta™** will remove hair conveniently and exfoliate skin utilizing specially designed exfoliating pads.

We worked long and hard to make this device suitable for your skin care and hair removal needs. After extensive testing, we added this product to our health & beauty collection for your home.

Enjoy your **Derma Seta™** ... and, once again, thank you!

**CONTENTS**

WARNING ..... 4 - 5

SAFETY PRECAUTIONS ..... 5 - 6

INTRODUCTION ..... 7

BENEFITS ..... 8

PARTS LIST - Basic Unit ..... 9

PARTS LIST - Deluxe Unit ..... 10

OPERATION ..... 11 - 14

*Assembly* ..... 11

        • *Basic* ..... 11 - 12

        • *Deluxe* ..... 13 - 14

TIPS ..... 15

ATTACHMENTS ..... 16 - 20

    • *Included Attachments* ..... 16

    • *Bonus Items* ..... 17

    • *Optional Items* ..... 18

    • *Changing Attachments* ..... 19


    • *Replacing Pads* ..... 20

TROUBLESHOOTING ..... 21

• **REPLACING DERMA CRYSTAL HAIR REMOVAL PADS**




Ensure **Derma Seta™** is turned OFF before replacing Derma Crystal Hair Removal Pads.

- 1) Peel the existing Derma Crystal Hair Removal Pad off the rotating head.
- 2) Take the new Derma Crystal Hair Removal Pad and peel off the protective sticker on the back.
- 3) Align the Derma Crystal Hair Removal Pad edge with the rotating pad's edge.

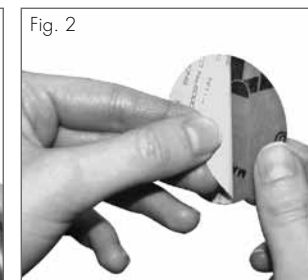

- **CHANGING ATTACHMENTS**

- 1) Ensure **Derma Seta™** is turned OFF before changing the attachment.
- 2) Remove by putting your finger in a small indentation on the side of the rotating pad. Then pull the rotating pad while your finger is placed in the indentation.
- 3) Pick desired attachment. Line up the slots on the head of the **Derma Seta™** device with the appropriate clips on the back of the attachment head. Gently push the attachment until secure.

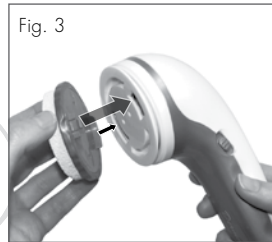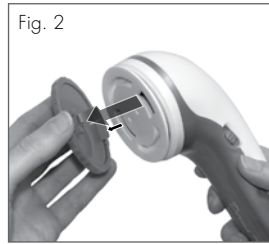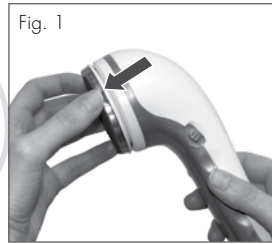

- **WARNING**

- Use this appliance only for its intended use as described in this manual.
- Do not use attachments not recommended by the manufacturer.
- Close supervision is necessary when this appliance is used by individuals with disabilities.
- Adult supervision necessary when used on or by children.
- Not for use on animals.
- Do not place appliance on any surface while it is operating.
- **Do not apply excessive pressure when performing the treatment. Applying too much pressure will reduce the effectiveness of the treatment. The motor is geared to either stop or slow significantly if too much pressure is applied.**
- If you are removing a lot of hair, the first treatment may take some time.

- **CAUTION**

- If you have a preexisting skin condition, sensitive skin or circulation problem, you should consult a dermatologist before using this device.
- For external use only.
- Do not remove hair from moles or areas of sunburn or abrasions.
- Avoid treating sunburned skin.
- Only use specified attachments for your bikini, neck and facial area. For more details please refer to the operation section of the manual.
- Do not apply lotion before use.
- Do not use if you suffer from skin conditions or irritations such as psoriasis, eczema, sensitive skin, skin allergies, rashes or circulatory problems. If you are at all unsure, consult a medical professional before use.

### **⚠ CAUTION** - continued

- Do not try to remove hair from open wounds, cuts, burns or warts.
- Do not use in the same area for a long period of time as it may cause a skin burn.

**! Note:** Carry out a patch test on a small area of your skin 24 hours before use. If irritation occurs discontinue use and consult your dermatologist.

**! Note:** Some people may experience a slight redness after use. This is quite normal and will fade within an hour or two. Ensure that your skin is completely dry before treatment.

### **SAFETY PRECAUTIONS**

As with any electrical device, please take the appropriate precautions when operating this beauty aid and read all instructions before using.

- 1) To reduce the risk of death or serious injury from electrical shock:
  - Do not immerse the device in water. Unplug it immediately if it comes in contact with water.
  - Do not use in the bath or shower.
  - Always unplug the device after using or before cleaning.
- 2) Use this appliance only for its intended use as described above and do not use any attachments not recommended by the manufacturer.

### • *OPTIONAL ITEMS*

#### • *Basic Unit Adaptor*

Once you are out of battery power, you can always plug in the adaptor for continuous use.

**! Caution:** Before using the adaptor please remove the batteries from the Basic Unit.



#### • *Deluxe Unit*

The unit also comes with rechargeable built in batteries and a charger, so you never have to worry about buying new batteries again.

The speed switch indicated in the part list lets you choose between 2 speeds, for beginner and advanced users.



- BONUS ITEMS

- Attachments

#### Foot Smoother

This coarse grit attachment is specifically designed for rough spots on your feet.

**! Note:** DO NOT use for hair removal; ONLY use for its intended purpose.

#### Macro & Micro Massage Modules

Choose from two sizes to help reduce cellulite. It also improves circulation for a smoother appearance.

#### Deep Cleanser Brush

Use to enhance the benefits of all face cleansers.



#### SAFETY PRECAUTIONS - continued

- 3) Never leave the device unattended when plugged in and keep the device away from any heated surfaces.
- 4) Never operate or recharge if the plug or electrical cord connected to the Charging Adaptor is damaged or has come in contact with water.
- 5) Do not drop or insert any foreign objects into the rotating head.
- 6) Never mix batteries. Replace all cells when weak.
- 7) Do not try to recharge non rechargeable batteries. Charging non rechargeable batteries may generate hydrogen that can lead to an explosion.
- 8) Observe the right polarity when installing or replacing batteries. Wrongly installed (reversed) battery can heat up and create a fire hazard.
- 9) Remove fully discharged batteries from the unit. Old cells can leak and cause corrosion.
- 10) Remove the batteries when the unit is not used for a longer period to prevent corrosion.
- 11) Always keep batteries away from children.

## INTRODUCTION

Welcome to the ultimate process in painless and easy hair removal. Now you can remove unwanted hair from arms and legs quickly and effectively, leaving your skin feeling smooth, silky and radiantly beautiful. With a variety of attachments, **Derma Seta™** provides the total skin-pampering experience in the comfort of your home or wherever you travel.

### • INCLUDED ATTACHMENTS

**Derma Seta™** comes equipped with a variety of convenient attachments to provide a total skin beauty experience:

**Derma Crystal Hair Removal Pad** – this fine grit attachment is used to remove hair and exfoliate skin at the same time. 5 additional exfoliating pads are included with the unit. Please see the instructions on how to replace Derma Crystal Hair Removal Pads on page 20.

**! Note:** ONLY use this attachment for hair removal and exfoliation.

**Pumice Stone** – attach this to remove calluses or rough skin on the feet. The Pumice Stone is to be used on wet, soaked feet for best results.

**Body Buffer** – attach to effectively apply lotion to your body as well as to remove dry skin after exfoliation.

**Facial Hair Remover** – use to effectively remove any unwanted facial hair.

**Bikini Line Hair Remover** – use to remove unnecessary hair in the bikini line area.





## TIPS

- Give yourself some time to determine the best motion for maximum results.
- Results may vary, but typically using **Derma Seta™** 10 to 12 times will achieve the optimum effect.
- The Derma Crystal Hair Removal Pad is intended for dry use only and will be damaged if used with any liquids, lotions, shaving creams or gels.
- When the pad becomes noticeably smoother, change the pad to achieve the best results.

## BENEFITS

We understand how hectic life can get and keeping up appearances can be very time consuming. **Derma Seta™** was developed with your busy schedule in mind. Our design is very user-friendly and extremely easy to operate, leaving you more time in the day to get to the most important things in life. Compared to shaving, waxing and pedicures, **Derma Seta™** takes less time to acquire the same results. By eliminating any potential for pain, additional costs, and time to visit a professional, **Derma Seta™** is the ideal product to get the results you want, when it's convenient for you.

**Derma Seta™** is designed to remove hair conveniently and exfoliate skin utilizing specially designed Derma Crystal Hair Removal Pads. These pads are excellent for gentle removal of hair on arms and legs and, unlike laser hair removal, is suitable for all skin types.

This multipurpose beauty aid is compact, portable and easy to use. With a variety of accessories, it is perfectly designed to remove calluses, smooth skin and apply lotions as well as routine hair removal.

Never again will you have to worry about nicks and abrasions from shaving, and **Derma Seta™** gives you effective and painless removal of unwanted hair. Leaving soft hairless skin and heightened blood circulation for 1-2 days, **Derma Seta™** is the definitive product for easy, pain-free hair removal.

PART LIST - BASIC UNIT



- 6) The attachment head will rotate clockwise only. Using light pressure (keeping the pads in one place for too long may cause skin injury) apply the exfoliating disc to your skin.
- 7) Move the rotating disc in a gentle circular motion in a counterclockwise direction, gradually covering the entire area. You will feel a warm sensation, which is common. With practice, you will determine the exact amount of pressure to effectively remove hair.
- 8) The exfoliating action produces a fine white powder accumulating on your skin. This is normal as the hair removal process reduces hair and dry skin to a powder that can be removed with the soft brush attachment. (See attachment section for further instructions).
- 9) Feel your skin as you move the pads. If you miss a spot, or encounter a stubborn area, gently pass over the area again in a circular motion, until the hair is gone.
  - ! **Caution:** Using Derma Seta™ for too long on one area may cause a skin burn.
- 10) Once finished, push the power switch to the OFF position and place back into the Charging Cradle.
  - ! **Note:** Make sure the unit is OFF before placing it in the stand.

• DELUXE UNIT

1) Attach the Charging Cradle to the adaptor provided. Plug the cord into an electrical outlet; a light will appear on the front of the Cradle. If the light is green, **Derma Seta™** is charged. If the light is red the unit still requires charging (for best results, **Derma Seta™** should be fully charged).

**! Note:** Charging the unit can take up to 4 hours.

2) Once **Derma Seta™** is charged (indicated by the green light on the Cradle) ensure power switch is in OFF position and remove from Cradle.

3) Peel backing from Derma Crystal Hair Removal Pads and affix the accessory onto the rotating head.

4) Line up the slots on the rotating head of the **Derma Seta™** device with the appropriate clips on the back of the attachment head. Gently push the attachment until secure.

5) Once attachment head is securely in place, push the power switch to the ON position (the rotating head will begin to spin). By pushing the power switch up one and two notches, you will increase the speed of the rotating head. You can pick your desired speed; Low and High.

**! Caution:** Do not apply lotion before use.

PART LIST - DELUXE UNIT



## OPERATION

Derma Seta™ and its attachments have all been carefully crafted to offer ease of operation coupled with excellent results. The grain on the pads has been specifically formulated to remove arm and leg hair.

*For the best results, please follow these steps:*

**Before use:** If you have a long hair style, make sure to keep it away from the rotating pads as it may become entangled. Do not use this device if you have any skin irritations or cuts.

## ASSEMBLY

### • BASIC UNIT

- 1) Open Battery Pack cover located on handle of Derma Seta™, and insert six alkaline AA Batteries (follow diagram inside Battery Pack). Replace cover.
- 2) Peel backing from Derma Crystal Hair Removal Pads and affix the accessory onto the rotating head.
- 3) Line up the slots on the head of the Derma Seta™ device with the appropriate clips on the back of the attachment until secure.

- 4) Once attachment head is securely in place, push the power switch (located on the side of the handle of the Derma Seta™ device) to the ON position.
- 5) The attachment head will rotate clockwise only. Using light pressure (keeping the pads in one place for too long may cause skin injury) apply the exfoliating disc to your skin.
- 6) Move the rotating disc in a gentle circular motion in a counterclockwise direction, gradually covering the entire area. You will feel a warm sensation, which is common. With practice, you will determine the exact amount of pressure to effectively remove hair.
- 7) The exfoliating action produces a fine white powder accumulating on your skin. This is normal as the hair removal process reduces hair and dry skin to a powder that can be removed with the soft brush attachment. (See attachment section for further instructions).
- 8) Feel your skin as you move the pads. If you miss a spot or encounter a stubborn area, gently pass over the area again in a circular motion until the hair is gone.

**! Caution:** Using Derma Seta™ for too long on one area may cause a skin burn.

- 9) Once finished, push the power switch to the OFF position.

**! Note:** Make sure the unit is OFF before placing it in the stand.

*DermaSeta™*

Distributed By:

Thane International, Inc.,  
La Quinta, CA 92253  
[www.thane.com](http://www.thane.com)

Thane Direct Canada Inc.,  
Mississauga, ON  
[www.thane.ca](http://www.thane.ca)

Thane Direct UK Ltd.,  
Admail ADM3996,  
London, W1T 1ZU  
[www.thanedirect.co.uk](http://www.thanedirect.co.uk)



Copyright © 2010 Thane International, Inc.

DS-REV 1-ENG 09/04/10